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
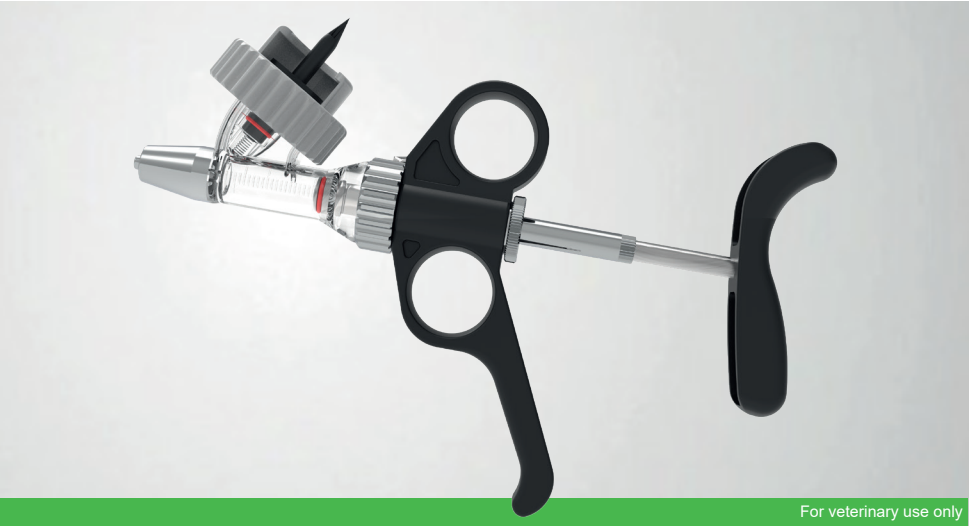
ENG

Automatic self-filling syringe

- 1) Ergonomic, rotating handle
- 2) Cylinder with bottle attach
- 3) Piston and O-ring
- 4) Dosage setting
- 5) Universal adapter
- 6) Tube adapter
- 7) Lock nut

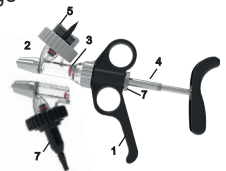The HSW B-MATIC is designed for robust use. Thanks to the special design of the handle, the syringe can be easily converted from a bottle attachment to a tube attachment, depending on the version. It is very easy to dismantle and clean.

Fitting the needle:

! For your safety, only remove the protective cap from the needle, when using disposable needles, immediately before using the syringe.

- Place the needle onto the Luer Lock fitting and secure it by turning the needle fitting clockwise..
- To remove the needle, follow the procedure in reverse. Refit the protective cap.

Connecting the tube:

- Loosen the lock nut so that the handle can be pushed backwards and the pin is released.
- Turn the cylinder to the required position.
- Push the handle forwards so that the pin secures the position, then tighten the lock nut.
- Attach the tube to the tube adapter and push it firmly into place.
- Connect the other end of the tube to the ascending tube or withdrawal cap.

Attaching the bottle:

- Select the appropriate universal adapter for the bottle: White: bottle neck diameter 20 mm, for 50–100 ml bottles, Blue: bottle neck diameter 30 mm, for 100–250 ml bottles, Green: bottle neck diameter 33 mm, for bottles larger than 500 ml

- Insert the bottle into the universal adapter.
- Insert the suction tip into the centre of the bottle.
- Tighten the universal adapter without turning the bottle.

Setting the dosage volume:

- Turn the lock nut anti-clockwise and loosen it.
- Turn the dosage sleeve clockwise or anti-clockwise until the required dose indication on the cylinder is reached. The front edge of the piston serves as the dosing edge.
- Tighten the lock nut.

Preparing for use:

- Remove the protective cap from the cannula.
- Hold the B-MATIC at a slight incline and make sure that the suction tip remains below the liquid level to prevent air from being drawn into the cylinder.
- Press the movable handle and allow it to return to draw liquid into the cylinder.
- If air bubbles are present, repeat this procedure until no air bubbles are visible.
- When filling the syringe and preparing the injection, make sure that the syringe is set to 0.5 ml for 0.5 ml syringes and to at least 0.5 ml for 1 ml syringes. Also ensure that there are no air bubbles in the cylinder.
- If necessary, point the needle upwards so that the air can escape. After this initial setting, any required dosage volume may be selected.

Maintenance after use:

- After application, remove any residue using standard cleaning agents and rinse the syringe thoroughly under warm water.

Disinfection:

- ! *It is strongly recommended to disinfect the product and needles before and after each use.*
- Clean the syringe with hot water.
- Place the entire syringe, including the attached needles, in a container of water and boil it for a maximum of 20 minutes.
- After boiling, remove the syringe, dry it with a cleaning cloth, and flush it with water until dry.

Lubrication:

- Lubricate the O-ring with medical white oil or, alternatively, with 100% pure extra virgin olive oil, cold-pressed, or an equivalent edible oil.
- Observe the manufacturer's shelf-life and storage instructions.
- Do not use mineral oils, silicone oils, WD-40, glycerine or alcohol. Petroleum-based oils may damage the piston or reduce the service life of the product.
- After lubrication, the B-MATIC can be reassembled by screwing the cylinder onto the handle. Loosen the adjusting screw on the piston rod and remove the cylinder.

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General notes and warnings for the HSW B-MATIC

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General warnings for all products of the HSW B-MATIC family:

- ! The HSW B-MATIC is only intended to be used for the treatment of animals. The application on humans may cause damage to their health.
- ! Do not point the HSW B-MATIC towards any part of your body or towards other people! This can lead to injuries or accidental injections.
- ! In case of accidental injection or application of a substance, consult a doctor immediately and present the information sheet of the respective substance to the doctor.
- ! The use of non-approved substances or the overdosage / underdosage of these approved substances can lead to adverse health effects in animals and/or later in humans.
- ! Caution, parts of the B-MATIC can break / be damaged. Damaged parts can cause injuries to the operator and animals.
- ! Caution, sharp needles and broken needles can cause harm to the operator and animals. Use needles with care and do not bend them.

General notes for all products of the HSW B-MATIC family:

- Substances to be administered with the HSW B-MATIC are not provided by HSW. The manufacturer's instructions must be adhered.
- The HSW B-MATIC may only be used by operators who have read and understood the user manual in full.
- It is recommended to wear protective goggles, a mouth guard and rubber gloves during the use and cleaning of the HSW B-MATIC.
- Transport and store in dry and ambient conditions.

Explanation of packaging symbols:

eIFU:

- This symbol indicates that a digital instruction manual is available for this syringe. The digital instruction manual can be found under the following link: www.henkesasswolf.de/manuals

Recycling codes:

- The recycling codes provide information about the material used and can be helpful for waste separation.
- The letter abbreviation PAP identifies cardboard and paper. PAP 20 = Corrugated cardboard

Disposal:

- The HSW B-MATIC products can be disposed of in household waste.
- Due to the risk of injury, cannulas / needles must be disposed of as hazardous waste in accordance with the regulations of each country.

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Intended Use & Intended User:

Intended Use: Intramuscular / subcutaneous vaccination & intranasal application of livestock

Intended User: Vaccination as well as maintenance and repair can be carried out by a veterinarian or instructed / trained farm staff.

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You can download an updated version of the complete spare part list for this product under www.henkesasswolf.de or kindly request the same.

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© copyright Henke-Sass, Wolf GmbH